

ergoflīx

Mobility for Your Everyday Life

ergoflīx[®] basic



ergoflīx

UP TO 100 %
REIMBURSEMENT



Contents

04 Flexibility

Practical folding mechanism

06 Comfort for Your Everyday Life

Comfortable from the very first moment

08 Intuitive Controls

For independent mobility

10 Energy for Your Day

Power that moves you

11 Your Safe Driving Experience

Reliable mobility wherever you go

12 The ergoflix® basic at a Glance

Key Features

15 Range of Accessories

For greater individuality

16 Health Insurance Funding

Save costs now

18 Book Your Test Drive

Free of charge and with
no obligation




Flexibility

Practical folding mechanism

The ergoflix® basic can be folded in seconds in just a few simple steps. This allows you to store it compactly and transport it in your car.

- Operate the lever on the back of the wheelchair and gently push the backrest forward.
- Grip the footrest and pull it towards the backrest. Your wheelchair is now compactly folded.

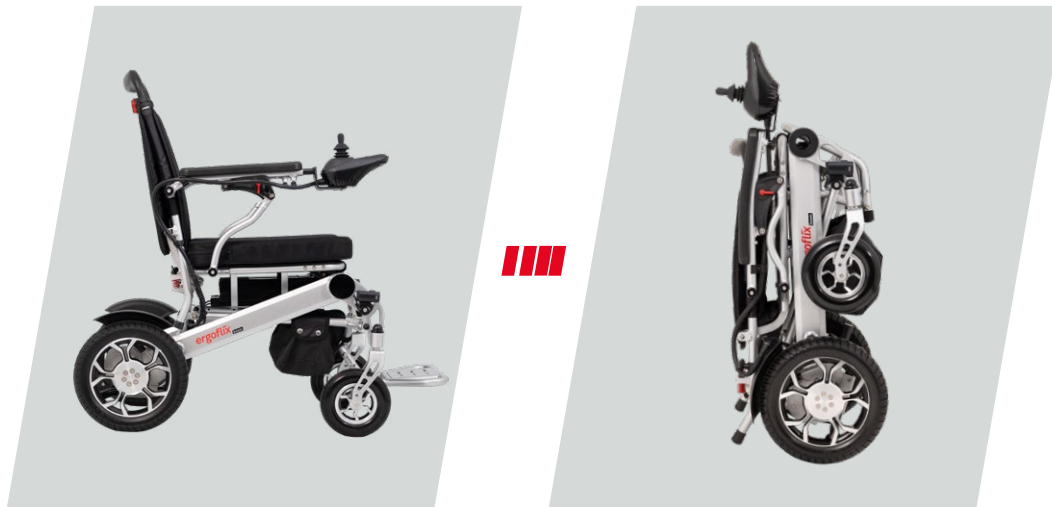



Comfort for Your Everyday Life

Comfortable from the very first moment

Flip-up armrests

- I** for easy access when getting in and out
- I** for moving close to tables or bathroom facilities



Adjustable firmness

- I** adjust the firmness of the backrest and seat
- I** three tension straps allow flexible adjustment to your body shape



Suspended seat

- /** reduces shocks and vibrations
- /** provides a comfortable ride even on uneven surfaces

8-inch front wheels

- /** puncture-proof tyres
- /** provide excellent manoeuvrability, especially indoors



Intuitive Controls

For independent mobility

Simply sit down and set off! The control module of the ergoflix® basic is operated intuitively using the joystick, allowing you to master the controls in no time. Enjoy a reassuring sense of safety from the very first moment.

Additional Features:

- Battery charge indicator
- Speed adjustment
- Conveniently charge your wheelchair via the control module





Energy for Your Day

Power that moves you

Motor with replaceable brushes

- low-maintenance
- reliable in everyday use

Maintenance-free lithium-ion battery

- for a range of up to 15 km¹ on a single charge
- easy to remove and charge separately



The range can be extended to up to 30 km¹ with a second battery.



Your Safe Driving Experience

So you can always reach your destination reliably

For added safety, the ergoflix® basic is equipped with additional features. You can rely on your mobile companion at all times.

■ The **lights at the front and rear** of the wheelchair provide good visibility and make you more visible in road traffic.

StVZO ✓



■ **Reflective elements** further enhance your visibility.

■ The **integrated anti-tip wheels** provide additional stability when travelling uphill and downhill.

The ergoflix[®] basic

Key Features at a Glance

**RANGE EXTENDABLE
TO UP TO 30KM¹**

Including lighting for
**optimum visibility
at any time of day**

Raised seating
position **for easy
access when
getting in and out**

12.5-inch rear wheels
are **puncture-proof and
easily overcome smaller
obstacles**

Suspended seat
absorbs shocks

Extendable range
with an additional
battery

The image shows the ergoflix[®]
basic with two batteries – one
battery is included as standard.

Low-maintenance brush
motors for **high reliability
in everyday use**



Performance and Range

Maximum range ¹	15 km (extendable to up to 30 km with an additional battery)
Maximum load capacity	120 kg
Motor	2 × low-maintenance 200 W brush motors
Gradient capability	14 %
Batteries & charging time (per battery)	1 × 24 V / 12 Ah approx. 6 h

Dimensions and Transport

Weight	22.7 kg 24.5 kg with battery
Dimensions (LxWxH)	103x61x99 cm
Folded dimensions (LxWxH)	38.5x61x79.5 cm
Seat (LxWxH)	43.5x48x55 cm

Comfort and Features

Suspension	Suspended seat
Puncture-proof tyres	
Front	PU 8"
Rear	PU 12.5"
Lighting	Standard lighting

**Health Insurance
Funding and Price**

Health insurance reimbursement code	18.50.07.0013
Price	€3,890

¹ Actual range may vary in practice depending on weight, driving style, speed, use of comfort and auxiliary functions, ambient temperature, payload and terrain.

Range of Accessories

For Greater Individuality

Telescopic ramp



Wheelchair bag | rear



Organiser bag



14

Protective cover with transport bag



FURTHER ACCESSORIES ARE
AVAILABLE ON REQUEST




Save up to 100% of the Cost

Secure Your Health Insurance Funding

1 Doctor's Prescription

The first step towards receiving health insurance funding is an appointment with your doctor. If your doctor confirms that you need a mobility aid, they will issue a medical prescription, similar to a prescription for medication. The prescription should include the diagnosis, e.g. "immobility" or "gait disorder".

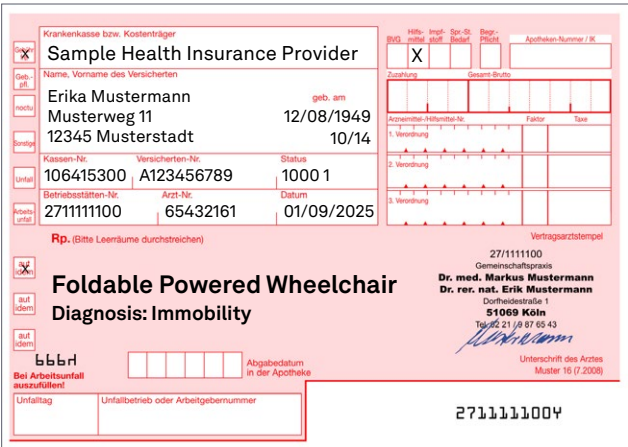

The image shows a sample German medical prescription form for a 'Foldable Powered Wheelchair'. The form is divided into several sections. The top left section contains patient information: Name (Erika Mustermann), birth date (12/08/1949), and address (Musterweg 11, 12345 Musterstadt). The top right section contains insurance information: Krankenkasse (Sample Health Insurance Provider), Kassen-Nr. (106415300), and Versicherten-Nr. (A123456789). The middle section contains the doctor's information: Dr. med. Markus Mustermann, Dr. rer. nat. Erik Mustermann, and their address (Dorfriedestraße 1, 51069 Köln). The bottom section contains the diagnosis (Immobility) and the date of the prescription (01/09/2025). The form also includes a section for the pharmacy (Apotheken-Nummer / IK) and a section for the patient's signature and date.

Sample Prescription

2 Contact ergoflix

Once you have received your prescription, please get in touch with us. Together, we will arrange a personal consultation at a mobility retailer near you. There, you can find the right ergoflix wheelchair for you.

3 Submitting Your Documents

The mobility retailer submits all documents to your health insurance provider. Your provider will review the application and then contact you directly. Once approval has been granted, you will receive your ergoflix and can contact your mobility retailer directly with any questions regarding your wheelchair.

Good to know: For Germany, the health insurance reimbursement code for the ergoflix® basic is 18.50.07.0013.

Note for Private Health Insurance Policyholders

The process for private health insurance differs from the statutory procedure. Here too, everything begins with a medical prescription and a consultation at a mobility retailer. You will then receive a quotation. Submit this quotation to your private health insurance provider. Once approval has been granted, you can order your ergoflix directly through your mobility retailer.

Do you have questions about your medical aid application, or would you like to appeal following a rejection? Send us an email at **krankenkasse@ergoflix.de** or simply call us on **+49 2852 9459094**.

We are happy to support you throughout the application process. Our experts are highly experienced in dealing with health insurance procedures.



Book Your Test Drive

Free of charge and with no obligation

Get to Know the ergoflix® basic


At ergoflix, you can test drive your future companion free of charge. At a time that suits you, one of our certified medical device consultants will visit you at home, explain the wheelchair's features and answer all your questions. You can then try out the wheelchair for yourself in your own home.

Alternatively, you are welcome to arrange a consultation at a mobility retailer near you. We work with more than 800 partners throughout Germany.

Booking Your Test Drive

Together, we will find a suitable appointment.

To book your test drive, simply call us on **+49 2852 9459000**.

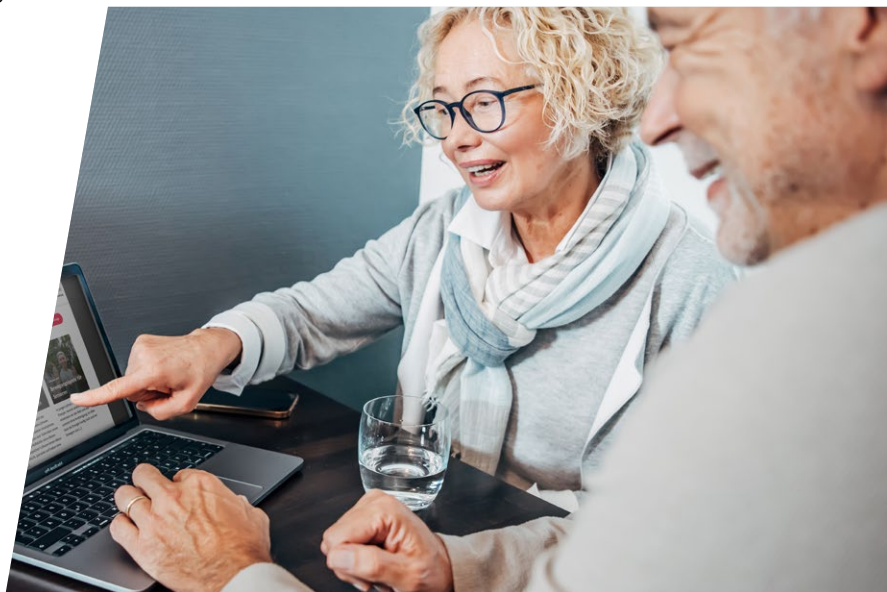



ergoflix Group GmbH
Marktstr. 20
46499 Hamminkeln
Phone +49 2852 9459000
Fax +49 2852 9459099
info@ergoflix.de
www.ergoflix.de

A Sunrise Medical® Brand  **SUNRISE**
MEDICAL

Visit us on social media



Telephone service hours

March to November:

Monday to Friday from 8:30 a.m. to 6:00 p.m.

December to February:

Monday to Friday from 8:30 a.m. to 5:00 p.m.

Manufacturer

Anhui JBH Medical Apparatus Co., Ltd.
No. 116 Qicang Road, Mingguang, Chuzhou
Anhui Province, Anhui, 239400
CHINA



SUNGO Europe B.V.
Fascinato Boulevard 522, Unit 1.7, 2909VA
Capelle aan den IJssel, The Netherlands

Importer and Distributor

ergoflix Group GmbH
Marktstr. 20
46499 Hamminkeln
GERMANY

PHONE +49 2852 9459000
FAX +49 2852 9459099

info@ergoflix.de
www.ergoflix.de

Legal notice: ergoflix Group GmbH acts as importer in accordance with Regulation (EU) 2017/745 (MDR). No liability accepted for printing errors. Errors and omissions excepted. Any colour deviations are due to the printing process. All dimensions are approximate values. All items shown without decoration. ergoflix Group GmbH acts as importer in accordance with the Medical Devices Act.

All prices are in EUR, including VAT and shipping within Germany.

Item no.: WM140



202608-06